



Activity Chair

In a few short years the Activity Chair has become the classroom standard for children and young adults with special positioning needs. Not only that, but parents who see what these chairs can do in school are getting them for use at home. We've also seen a growing number of speech and feeding specialists adopt this chair in their hospital-based clinics. It's safe to say the Activity Chair has truly revolutionized adaptive seating and continues to delight therapists with its versatility for people of all shapes and sizes, all conditions and special needs.



Find inservice videos, letters of medical necessity and informative articles at:
www.rifton.com/activitychair

"These chairs are the talk of the school. I've been trying as many kids as I can each day, but no one wants to give the chair back! The adjustments are fabulous, but it's how comfortable, functional, and straight the kids are in the chair that really excites us. *"*

*Judi Rogers, OTR/L
LAUSD, California*

● **Hi/lo Activity Chair (R870)** with
contoured headrest, laterals and tray.

DESIGNED FOR USE WITH THE
MOVE
PROGRAM

Activity Chair features



Spring option
The springs on the standard base (in both backrest and seat) allow for self-generated, dynamic motion—ideal for clients with autism. The spring column can be locked when not in use.

Standard base



Caster swivel lock
The swivel lock prevents the chair from drifting sideways when it is pushed for longer distances.



Adjustable back



Tilt-in-space



Hi/lo base



Adjustable back



Tilt-in-space



Adjust while in use

The Activity Chair can be easily adjusted without tools, with a client in the chair. This saves time and makes adjustments much simpler to gauge.



Easy height adjustment

The large and medium hi/lo bases are height adjustable with the foot pump.

Footboard features



Knee angle



Footplate angle



Flip up footboard

Perfect for the MOVE™ program



When used with forearm prompts, the Activity Chair becomes a forward-leaning instructional seating device as used in the MOVE™ program.

Transfers made easy



Remove the armrests and raise the hi/lo chair to the ideal transfer height for secure, safe lateral moves.



Flip the footboard up and tilt the chair forward for easy sit-to-stand transfers.



Lower the hi/lo chair until the footboard is flat on the floor. Because the footboard is strong enough for stand and pivot transfers, clients can assist in their own transfers.

A chair for everything

The Activity Chair's tilt-in-space feature allows easy transitions between forward leaning for task participation and unimpeded interaction, and reclining for rest or a higher level of support.



All Rifton Activity Chairs have tilt-in-space.

Whether you choose the standard or hi/lo base, you'll get Rifton's quick-adjusting, tilt-in-space feature. The standard base tilts 15° forward and 15° back, while the hi/lo tilts 15° forward and 25° back. Tilt-in-space is a valuable feature because it facilitates sit-to-stand and offers a variety of resting positions. Best of all, our tilt-in-space can be activated while your client is in the chair.





This medium standard chair is in posterior tilt to provide support for a student with poor head control.



The tilt-in-space feature accommodates users who need a higher degree of head support.

Infinitely adjustable for every need

The favorite of feeding clinics

Rifton's Activity Chair can be either a basic chair for a client who needs just a little support, or it can offer significant support for more involved individuals. The hi/lo base allows the client to be fed in a quiet, individual setting or at the family dinner table. Because you can measure the angle of the tilt-in-space and backrest recline, home caregivers can exactly duplicate the angles recommended by therapists.

Each size of Activity Chair has so much growth potential that it will provide many years of service. And it is easy to clean—essential for a feeding chair. The tray clips on and off with one hand.

Quickly adjustable tilt-in-space during feeding is beneficial. Therapists may start out the session doing oral motor exercise with the client upright and then tilt back to reduce the need for trunk control. They can then adjust the headrest forward to keep the head in a neutral position.

All these features make it the ideal feeding chair.



From an activity at counter height to a low-table group activity, the hi/lo base chair makes participation easy.

Components



Choice of arm supports

Armrests provide lateral boundaries as well as upper extremity weight-bearing assist to aid trunk control. Rifton's armrests can be raised or lowered individually and can be set at angles from horizontal to 30° up or down. The tray attaches to the armrests.

small 5.5 x 22.5 cm, **large** 8 x 30 cm



Forearm prompts facilitate head and trunk control while sitting, and secure distal placement of upper extremities. This enables the use of the pelvic girdle and trunk core muscles for improved postural stabilization and strengthening. Rifton's forearm prompts are fully adjustable for optimal positioning.

Padded area: small 10 x 18 cm, large 11 x 20 cm

Range of adjustment from elbow to fist: small 20–32 cm, large 27–41 cm



Choice of adjustable legs

Adjustable legs (long and short): Height-adjustable in 3 cm increments, these legs allow you to position your client at optimal height. The grippy tips of the stationary legs keep the chair from sliding and moving across the floor, even with dynamic rocking.

Adjustable legs with casters (long and short): Height-adjustable in 3 cm increments, these legs allow you to position your client at optimal height. The casters make your standard base chair easy to move when transitioning from one classroom to another.



Pads: Comfortable seat and backrest pads are essential for providing cushioned support and preventing pressure sores. The backrest and seat pads are made of polyurethane foam covered with a cleanable, healthcare-grade polyurethane-coated fabric. The seat pad can be turned front to back for longer wear. Pads are available in six colors.



Components (continued)



Headrests optimally position clients with poor head control. Choose a flat, contoured or winged headrest. The winged headrest adjusts from almost flat to 90°. All headrests adjust front to back and up and down using Rifton's unique two-knob adjustment system.

Flat and contoured: 23 x 13 cm

Winged: 41 x 13 cm when flat to 15 cm between wings at 90°



Laterals: Necessary for clients who need good upper body support, laterals are important for proper alignment of the spine. Laterals have independent vertical, lateral and angle adjustment.

Small 13 x 10 cm, **large** 17 x 13 cm



Laterals with chest strap: With the addition of a wide chest strap to the laterals, clients with limited upper torso control gain added trunk support and security.



Push handles: Recommended for all chairs with wheels, push handles make it easy for a caregiver to maneuver and move a client in the Activity Chair from one place to another.



Tray: The tray attaches to the armrests and provides a work, play and feeding surface with raised edges. It can be set at angles ranging from horizontal to 30° up or down, and can be raised or lowered to meet your client's needs.

Small 48 cm wide x 46 cm deep, **medium** 56 cm wide x 51 cm deep

Large 65 cm wide x 57 cm deep



Handhold: The handhold attaches to the tray for arm and hand positioning. Handholds are especially useful for clients with spastic muscle tone and are helpful for proper joint alignment.





Chest strap: The wide and comfortably padded chest strap provides added trunk support.

Max. circumference: small 86 cm, medium 102 cm, large 132 cm

Width of strap: small 8 cm, medium 9.5 cm, large 10 cm



Butterfly harness: The butterfly harness provides maximum anterior support giving security, safety and freedom of movement.

Slim-cut butterfly harness performs the same positioning function as the regular butterfly harness. The slim cut butterfly harness is an alternative especially designed for larger female clients.



Thigh belt: The thigh belt provides knee adduction and additional support and security for the client's thighs.

Max. circumference: small 81 cm, medium 104 cm, large 132 cm



Pelvic harness: An alternative to a typical seat belt, Rifton's unique pelvic harness firmly positions a client's pelvis by securing hips and upper thighs without pressure on the abdomen. The pelvic harness provides a stable, comfortable base for postural control and encourages proper alignment of the spine. It also controls extreme extensor spasticity that can result in sacral sitting.

Max. thigh circumference: small 56 cm, medium 64 cm, large 79 cm



Hip guides attach to the arm supports. Placed at the sides of the pelvis to align a client's torso, they provide maximum support of trunk control while sitting. Hip guides can also be important for proper alignment of the spine. They independently adjust vertically and horizontally. For transfers, hip guides remove with the arm supports.

Dimensions: small 8.5 x 14 cm, large 10 x 23.5 cm



Components (continued)



Abductor: The abductor separates the knees and helps to control spasticity that causes leg scissoring. The abductor can be adjusted forward and back in three positions.

Dimensions: small 8 x 15 x 9 cm, large 10 x 19 x 10 cm



Adductors: Adductors provide a comfortable lateral boundary for a client's knees to maintain proper body and joint alignment.

One size only: 13 cm high x 15 cm long



Leg prompts are helpful for clients with limited muscle control. They maintain body and joint alignment when sitting by supporting each leg individually with a strap that passes around the front. Leg prompts can be adjusted in and out and laterally.

Max. circumference: small 36 cm, medium 41 cm, large 46 cm



Footboard: The footboard facilitates weight-bearing through the legs and feet during functional sitting activities, and can be used with sandals or ankle straps to meet specific positioning needs. The footboard can flip up, swing back out of the way for sit-to-stand transfers, or can be used as a step for transfers.



Ankle straps (require footboard) are helpful for clients with limited muscle control and excursion movement. They allow a moderate range of movement while improving independence in postural control.

Note: Ankle straps cannot be used together with sandals.

Max. circumference: small 30 cm, large 42 cm



Sandals (require footboard): For clients who are unable to control leg movement, sandals adjust to position client's feet correctly for optimal body alignment.

Note: Sandals cannot be used together with Ankle straps.


Small 8 x 18 cm, **medium** 10 x 23 cm, **large** 13 x 27 cm





Wedges (require sandals) can accommodate clients with specific ankle positioning needs or uneven leg lengths. Use them to raise heels or toes, or combine two wedges under one sandal to add flat lift.

Thickness: small 3 cm, medium 3.5 cm, large 4 cm

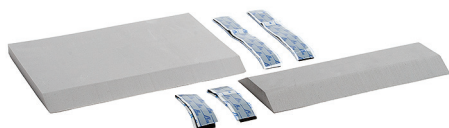

Lumbar and seat support kit: This kit lets you adapt the contours of the seat and back for a specific client's needs. The foam inserts can be cut to length and attached where needed under the seat and back pads with adhesive-backed hook and loop strips.



Backrest filler pad fills the space between the bottom edge of the backrest and the top edge of the seat when backrest is in one of the top three positions. It provides additional lower back support.



Mini kit: Consisting of a backrest insert pad and a mini trunk support, the mini kit supports children as short as 71 cm tall. It allows you to begin positioning children correctly at an earlier age, and continue to use the small chair for years before it's outgrown.

Dimensions: Trunk support max. circumference is 66 cm, backrest pad is 3 cm thick



Whitmyer adapter: Some Whitmyer headrests can be installed on all Rifton Activity Chairs. This adapter gives you the option to use a Whitmyer M2100-WB headrest mounting bracket and the Whitmyer Onyx headrest support system.

Note: Rifton supplies only the adapter and T20 power tip. Whitmyer mounting brackets and headrests must be purchased from a Whitmyer supplier.



Room to grow

Small chair with mini kit at its smallest dimensions

- 1 Backrest height 32 cm
- 2 Trunk support width 14 cm
- 3 Seat depth 18 cm
- 4 Seat width 18 cm
- 5 Armrest height above seat 13 cm
- 6 Seat height above footboard 18 cm




Small chair at its largest dimensions

- 1 Backrest height 39 cm
- 2 Distance between laterals 29 cm
- 3 Seat depth 30 cm
- 4 Seat width 30 cm
- 5 Armrest height above seat 19 cm
- 6 Seat height above footboard 30 cm




Small standard base
R820 Rifton Activity Chair

User dimensions (cm)	
Height	81–122
with mini kit	71–97
 Key user dimension: height The user's overall height is a general guide to help you select the appropriate chair. Choose the model that allows for growth. Important: Make sure that seat width, depth and height are adequate for each individual user, and that the user's weight does not exceed the maximum working load.	
Item dimensions (cm)	R820 small Standard base
Frame width	short legs: 53 long legs: 58
	short legs w/ casters: 57 long legs w/ casters: 60
Seat height above floor	short legs: 24–32 long legs: 47–55
	short legs w/ casters: 34–42 long legs w/ casters: 47–55
Seat angle (tilt-in-space)–degrees	15° forward, 15° back
Backrest angle–degrees	5° forward, 20° back
Footboard knee angle–degrees	45° – 110°
Seat height above footboard	18–30
Seat width with hip guides without hip guides	18–23 30
Seat depth with mini kit	20–30 18–28
Armrest height above seat	13–19
Trunk support width	14–29
Backrest height	32–39
Headrest height above seat	37–53
Max. working load (kg)	34



Small hi/lo base
R830 Rifton Activity Chair

User dimensions (cm)	
Height	81–122
with mini kit	71–97
 Key user dimension: height The user's overall height is a general guide to help you select the appropriate chair. Choose the model that allows for growth. Important: Make sure that seat width, depth and height are adequate for each individual user, and that the user's weight does not exceed the maximum working load.	
Item dimensions (cm)	R830 small Hi/lo base
Frame width	66
Seat height above floor	25–64
Seat angle (tilt-in-space)–degrees	15° forward, 25° back
Backrest angle–degrees	5° forward, 20° back
Footboard knee angle–degrees	45°–110°
Seat height above footboard	18–30
Seat width with hip guides without hip guides	18–23 30
Seat depth with mini kit	20–30 18–28
Armrest height above seat	13–19
Trunk support width	14–29
Backrest height	32–39
Headrest height above seat	37–53
Max. working load (kg)	34

2025 order form for small

Use dimension chart to select appropriate size.

		Standard R820 small	Hi/lo R830 small
* Seat and back includes seatbelt (Hi/lo includes footboard and push handles)	Adjustable	R821 <input type="checkbox"/>	R831 <input type="checkbox"/>
	Adjustable w/ back spring	R822 <input type="checkbox"/>	R832 <input type="checkbox"/>
* Base	Without spring	R823 <input type="checkbox"/>	R833 <input type="checkbox"/>
	With spring	R824 <input type="checkbox"/>	N/A
* Arm supports (pair)	Armrests	R827 <input type="checkbox"/>	R827 <input type="checkbox"/>
	Forearm prompts	R809 <input type="checkbox"/>	R809 <input type="checkbox"/>
* Adjustable legs	Short	R803 <input type="checkbox"/>	N/A
	Long	R894 <input type="checkbox"/>	N/A
	Short with casters	R805 <input type="checkbox"/>	N/A
	Long with casters	R896 <input type="checkbox"/>	N/A
* Pads	Tan	R815 <input type="checkbox"/>	R815 <input type="checkbox"/>
	Blue	R816 <input type="checkbox"/>	R816 <input type="checkbox"/>
	Red	R817 <input type="checkbox"/>	R817 <input type="checkbox"/>
	Pink	R971 <input type="checkbox"/>	R971 <input type="checkbox"/>
	Purple	R972 <input type="checkbox"/>	R972 <input type="checkbox"/>
	Green	R973 <input type="checkbox"/>	R973 <input type="checkbox"/>
Headrest	Contoured	R897 <input type="checkbox"/>	R897 <input type="checkbox"/>
	Flat	R898 <input type="checkbox"/>	R898 <input type="checkbox"/>
	Adjustable winged	R881 <input type="checkbox"/>	R881 <input type="checkbox"/>
Laterals (pair)	Small	R828 <input type="checkbox"/>	R828 <input type="checkbox"/>
	Small with chest strap	R819 <input type="checkbox"/>	R819 <input type="checkbox"/>
Push handles (recommended for all chairs with wheels)		R878 <input type="checkbox"/>	Included
Tray (requires armrests)		R667 <input type="checkbox"/>	R667 <input type="checkbox"/>
Handhold		R649 <input type="checkbox"/>	R649 <input type="checkbox"/>
Chest strap (wide)		R802 <input type="checkbox"/>	R802 <input type="checkbox"/>
Butterfly harness		R834 <input type="checkbox"/>	R834 <input type="checkbox"/>
Thigh belt		R811 <input type="checkbox"/>	R811 <input type="checkbox"/>
Pelvic harness		R845 <input type="checkbox"/>	R845 <input type="checkbox"/>
Hip guides (pair)		R814 <input type="checkbox"/>	R814 <input type="checkbox"/>
Abductor		R884 <input type="checkbox"/>	R884 <input type="checkbox"/>
Adductors (pair)		R892 <input type="checkbox"/>	R892 <input type="checkbox"/>
Leg prompt		R854 <input type="checkbox"/>	R854 <input type="checkbox"/>
Footboard (recommended for chairs with long legs)		R847 <input type="checkbox"/>	Included
Mini kit (backrest insert and trunk support)		R812 <input type="checkbox"/>	R812 <input type="checkbox"/>
Ankle straps (pair)		R899 <input type="checkbox"/>	R899 <input type="checkbox"/>
Sandals (pair)		R944 <input type="checkbox"/>	R944 <input type="checkbox"/>
Wedges (pair, require sandals)		R697 <input type="checkbox"/>	R697 <input type="checkbox"/>
Lumbar and seat support kit		R886 <input type="checkbox"/>	R886 <input type="checkbox"/>
Backrest filler pad		R849 <input type="checkbox"/>	R849 <input type="checkbox"/>
Whitmyer headrest adapter		R804 <input type="checkbox"/>	R804 <input type="checkbox"/>

* Required

Room to grow

Medium chair at its smallest dimensions

- 1 Backrest height **39 cm**
- 2 Distance between laterals **17 cm**
- 3 Seat depth **28 cm**
- 4 Seat width **22 cm**
- 5 Armrest height above seat **18 cm**
- 6 Seat height above footboard **19.5 cm**



Medium chair at its largest dimensions

- 1 Backrest height **48 cm**
- 2 Distance between laterals **30 cm**
- 3 Seat depth **41 cm**
- 4 Seat width **36 cm**
- 5 Armrest height above seat **27 cm**
- 6 Seat height above footboard **38.5 cm**



Medium standard base
R840 Rifton Activity Chair

User dimensions (cm)	
Height	102–142
	Key user dimension: height The user's overall height is a general guide to help you select the appropriate chair. Choose the model that allows for growth.
	Important: Make sure that seat width, depth and height are adequate for each individual user, and that the user's weight does not exceed the maximum working load.
Item dimensions (cm)	R840 medium Standard base
Frame width	short legs: 59
	long legs: 62
	short legs w/ casters: 62
	long legs w/ casters: 65
Seat height above floor	short legs: 32–42
	long legs: 48–58
	short legs w/ casters: 41–51
	long legs w/ casters: 48–58
Seat angle (tilt-in-space)–degrees	15° forward, 15° back
Backrest angle–degrees	5° forward, 20° back
Footboard knee angle–degrees	45°–110°
Seat height above footboard	19.5–38.5
Seat width with hip guides without hip guides	22–29
	36
Seat depth	28–41
Armrest height above seat	18–27
Distance between laterals	17–30
Backrest height	39–48
Headrest height above seat	44–62
Max. working load (kg)	68



Medium hi/lo base
R850 Rifton Activity Chair

User dimensions (cm)	
Height	102–142
	Key user dimension: height The user's overall height is a general guide to help you select the appropriate chair. Choose the model that allows for growth.
	Important: Make sure that seat width, depth and height are adequate for each individual user, and that the user's weight does not exceed the maximum working load.
Item dimensions (cm)	R850 medium Hi/lo base
Frame width	70
Seat height above floor	30–58
Seat angle (tilt-in-space)–degrees	15° forward, 25° back
Backrest angle–degrees	5° forward, 20° back
Footboard knee angle–degrees	45°–110°
Seat height above footboard	19.5–38.5
Seat width with hip guides without hip guides	22–29
	36
Seat depth	28–41
Armrest height above seat	18–27
Distance between laterals	17–30
Backrest height	39–48
Headrest height above seat	44–62
Max. working load (kg)	68

2025 order form for medium

Use dimension chart to select appropriate size.

		Standard R840 medium	Hi/lo R850 medium
* Seat and back includes seatbelt (Hi/lo includes footboard and push handles)	CHOOSE Adjustable	R841 <input type="checkbox"/>	R851 <input type="checkbox"/>
	Adjustable w/ back spring	R842 <input type="checkbox"/>	R852 <input type="checkbox"/>
* Base	CHOOSE Without spring	R843 <input type="checkbox"/>	R853 <input type="checkbox"/>
	With spring	R844 <input type="checkbox"/>	N/A
* Arm supports (pair)	CHOOSE Armrests	R868 <input type="checkbox"/>	R868 <input type="checkbox"/>
	Forearm prompts	R869 <input type="checkbox"/>	R869 <input type="checkbox"/>
* Adjustable legs	CHOOSE Short	R893 <input type="checkbox"/>	N/A
	Long	R894 <input type="checkbox"/>	N/A
	Short with casters	R895 <input type="checkbox"/>	N/A
	Long with casters	R896 <input type="checkbox"/>	N/A
* Pads	CHOOSE Tan	R835 <input type="checkbox"/>	R835 <input type="checkbox"/>
	Blue	R836 <input type="checkbox"/>	R836 <input type="checkbox"/>
	Red	R837 <input type="checkbox"/>	R837 <input type="checkbox"/>
	Pink	R974 <input type="checkbox"/>	R974 <input type="checkbox"/>
	Purple	R975 <input type="checkbox"/>	R975 <input type="checkbox"/>
	Green	R976 <input type="checkbox"/>	R976 <input type="checkbox"/>
Headrest	CHOOSE Contoured	R897 <input type="checkbox"/>	R897 <input type="checkbox"/>
	Flat	R898 <input type="checkbox"/>	R898 <input type="checkbox"/>
	Adjustable winged	R881 <input type="checkbox"/>	R881 <input type="checkbox"/>
Laterals (pair)	CHOOSE Small	R828 <input type="checkbox"/>	R828 <input type="checkbox"/>
	Small with chest strap	R829 <input type="checkbox"/>	R829 <input type="checkbox"/>
	Large	R838 <input type="checkbox"/>	R838 <input type="checkbox"/>
Push handles (recommended for all chairs with wheels)		R878 <input type="checkbox"/>	Included
Tray (requires armrests)		R668 <input type="checkbox"/>	R668 <input type="checkbox"/>
Handhold		R649 <input type="checkbox"/>	R649 <input type="checkbox"/>
Chest strap (wide)		R863 <input type="checkbox"/>	R863 <input type="checkbox"/>
Butterfly harness	CHOOSE Regular	R864 <input type="checkbox"/>	R864 <input type="checkbox"/>
	Slim cut	R882 <input type="checkbox"/>	R882 <input type="checkbox"/>
Thigh belt		R865 <input type="checkbox"/>	R865 <input type="checkbox"/>
Pelvic harness		R866 <input type="checkbox"/>	R866 <input type="checkbox"/>
Hip guides (pair)		R867 <input type="checkbox"/>	R867 <input type="checkbox"/>
Abductor	CHOOSE Small	R884 <input type="checkbox"/>	R884 <input type="checkbox"/>
	Large	R891 <input type="checkbox"/>	R891 <input type="checkbox"/>
Adductors (pair)		R892 <input type="checkbox"/>	R892 <input type="checkbox"/>
Leg prompt		R858 <input type="checkbox"/>	R858 <input type="checkbox"/>
Footboard (recommended for chairs with long legs)		R848 <input type="checkbox"/>	Included
Ankle straps (pair)		R899 <input type="checkbox"/>	R899 <input type="checkbox"/>
Sandals (pair)		R954 <input type="checkbox"/>	R954 <input type="checkbox"/>
Wedges (pair, require sandals)		R698 <input type="checkbox"/>	R698 <input type="checkbox"/>
Lumbar and seat support kit		R887 <input type="checkbox"/>	R887 <input type="checkbox"/>
Backrest filler pad		R849 <input type="checkbox"/>	R849 <input type="checkbox"/>
Whitmyer headrest adapter		R804 <input type="checkbox"/>	R804 <input type="checkbox"/>
Directional assist wheel retrofit kit		N/A	R965 (1 kit) <input type="checkbox"/> R966 (5 kits) <input type="checkbox"/>

* Required

Room to grow

Large chair at its smallest dimensions

- 1 Backrest height 48 cm
- 2 Distance between laterals 24 cm
- 3 Seat depth 38 cm
- 4 Seat width 28 cm
- 5 Armrest height above seat 18 cm
- 6 Seat height above footboard 28.5 cm



Large chair at its largest dimensions

- 1 Backrest height 61 cm
- 2 Distance between laterals 36 cm
- 3 Seat depth 51 cm
- 4 Seat width 43 cm
- 5 Armrest height above seat 27 cm
- 6 Seat height above footboard 54 cm




Large standard base
R860 Rifton Activity Chair



Large hi/lo base
R870 Rifton Activity Chair

User dimensions (cm)

Height 127–188

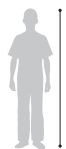

Key user dimension: height

The user's overall height is a general guide to help you select the appropriate chair. Choose the model that allows for growth.

Important: Make sure that seat width, depth and height are adequate for each individual user, and that the user's weight does not exceed the maximum working load.

Item dimensions (cm)

R860 large Standard base

Frame width	short legs: 68 long legs: 68
	short legs w/ casters: 70.5 long legs w/ casters: 70.5
Seat height above floor	short legs: 42–52 long legs: 51–61
	short legs w/ casters: 42–52 long legs w/ casters: 51–61
Seat angle (tilt-in-space)–degrees	15° forward, 15° back
Backrest angle–degrees	5° forward, 20° back
Footboard knee angle–degrees	45°–110°
Seat height above footboard	28.5–54
Seat width with hip guides without hip guides	28–36 43
Seat depth	38–51
Armrest height above seat	18–27
Distance between laterals	24–36
Backrest height	48–61
Headrest height above seat	50–75
Max. working load (kg)	113

User dimensions (cm)

Height 127–188



Key user dimension: height

The user's overall height is a general guide to help you select the appropriate chair. Choose the model that allows for growth.

Important: Make sure that seat width, depth and height are adequate for each individual user, and that the user's weight does not exceed the maximum working load.

Item dimensions (cm)

R870 large Hi/lo base

Frame width	75.5
Seat height above floor	41–61
Seat angle (tilt-in-space)–degrees	15° forward, 25° back
Backrest angle–degrees	5° forward, 20° back
Footboard knee angle–degrees	45°–110°
Seat height above footboard	28.5–54
Seat width with hip guides without hip guides	28–36 43
Seat depth	38–51
Armrest height above seat	18–27
Distance between laterals	24–36
Backrest height	48–61
Headrest height above seat	50–75
Max. working load (kg)	113

2025 order form for large

Use dimension chart to select appropriate size.

		Standard R860 large		Hi/lo R870 large			
* Seat and back includes seatbelt (Hi/lo includes footboard and push handles)	CHOOSE	Adjustable	R861	<input type="checkbox"/>	R871	<input type="checkbox"/>	
		Adjustable w/ back spring	R862	<input type="checkbox"/>	R872	<input type="checkbox"/>	
* Base	CHOOSE	Without spring	R873	<input type="checkbox"/>	R883	<input type="checkbox"/>	
		With spring	R874	<input type="checkbox"/>	N/A		
* Arm supports (pair)	CHOOSE	Armrests	R868	<input type="checkbox"/>	R868	<input type="checkbox"/>	
		Forearm prompts	R879	<input type="checkbox"/>	R879	<input type="checkbox"/>	
* Adjustable legs	CHOOSE	Short	R808	<input type="checkbox"/>	N/A		
		Long	R894	<input type="checkbox"/>	N/A		
		Short with casters	R895	<input type="checkbox"/>	N/A		
		Long with casters	R896	<input type="checkbox"/>	N/A		
* Pads	CHOOSE	Tan	R875	<input type="checkbox"/>	R875	<input type="checkbox"/>	
		Blue	R876	<input type="checkbox"/>	R876	<input type="checkbox"/>	
		Red	R877	<input type="checkbox"/>	R877	<input type="checkbox"/>	
		Pink	R977	<input type="checkbox"/>	R977	<input type="checkbox"/>	
		Purple	R978	<input type="checkbox"/>	R978	<input type="checkbox"/>	
		Green	R979	<input type="checkbox"/>	R979	<input type="checkbox"/>	
Headrest	CHOOSE	Contoured	R897	<input type="checkbox"/>	R897	<input type="checkbox"/>	
		Flat	R898	<input type="checkbox"/>	R898	<input type="checkbox"/>	
		Adjustable winged	R881	<input type="checkbox"/>	R881	<input type="checkbox"/>	
Laterals (pair)	CHOOSE	Small	R828	<input type="checkbox"/>	R828	<input type="checkbox"/>	
		Large	R838	<input type="checkbox"/>	R838	<input type="checkbox"/>	
		Large with chest strap	R839	<input type="checkbox"/>	R839	<input type="checkbox"/>	
Push handles (recommended for all chairs with wheels)		R878	<input type="checkbox"/>	Included			
Tray (requires armrests)		R669	<input type="checkbox"/>	R669		<input type="checkbox"/>	
Handhold		R649	<input type="checkbox"/>	R649		<input type="checkbox"/>	
Chest strap (wide)		R807	<input type="checkbox"/>	R807		<input type="checkbox"/>	
Butterfly harness	CHOOSE	Regular	R889	<input type="checkbox"/>	R889	<input type="checkbox"/>	
		Slim cut	R826	<input type="checkbox"/>	R826	<input type="checkbox"/>	
Thigh belt		R813	<input type="checkbox"/>	R813		<input type="checkbox"/>	
Pelvic harness		R846	<input type="checkbox"/>	R846		<input type="checkbox"/>	
Hip guides (pair)		R867	<input type="checkbox"/>	R867		<input type="checkbox"/>	
Abductor	CHOOSE	Small	R884	<input type="checkbox"/>	R884	<input type="checkbox"/>	
		Large	R891	<input type="checkbox"/>	R891	<input type="checkbox"/>	
Adductors (pair)		R892	<input type="checkbox"/>	R892		<input type="checkbox"/>	
Leg prompt		R859	<input type="checkbox"/>	R859		<input type="checkbox"/>	
Footboard (recommended for chairs with long legs)		R857	<input type="checkbox"/>	Included			
Ankle straps (pair)		R855	<input type="checkbox"/>	R855		<input type="checkbox"/>	
Sandals (pair)		R964	<input type="checkbox"/>	R964		<input type="checkbox"/>	
Wedges (pair, require sandals)		R699	<input type="checkbox"/>	R699		<input type="checkbox"/>	
Lumbar and seat support kit		R888	<input type="checkbox"/>	R888		<input type="checkbox"/>	
Backrest filler pad		R856	<input type="checkbox"/>	R856		<input type="checkbox"/>	
Whitmyer headrest adapter		R804	<input type="checkbox"/>	R804		<input type="checkbox"/>	
Directional assist wheel retrofit kit		N/A		R965 (1 kit)	<input type="checkbox"/>	R966 (5 kits)	<input type="checkbox"/>

* Required