



# E-Pacer

**The Rifton E-Pacer** combines a powerful electric lift column and trunk support with the key gait training features of the classic Pacer. Its sit-to-stand functionality removes a major barrier to gait training for larger or highly dependent clients. The optional thigh straps provide seat-to-seat transfer capability for added versatility.

- Zero lift sit-to-stand and seat-to-seat transfers
- Expandable base frame
- Accommodates users up to 6'5" and 350 lb



Scan here to order, request a quote, or see detailed product information.



# Key features

## 1 Zero lift transfers

The E-Pacer can perform sit-to-stand transfers as well as seat-to-seat transfers with no lifting on the part of the caregiver.



## 2 Expandable legs

Widen the base legs for clients needing extra step width or to get up close to a wheelchair prior to a transfer.



## 3 Height range

The 26" adjustment range of the E-Pacer's electric lift column provides versatility for use with clients of different statures.



## 4 Built-in scale

Optional Bluetooth scale and Gait Tracker app allow caregivers to track weight-bearing progress over time.





Scan for sizing and pricing.



# E-Pacer components



## Scale

Weigh clients during transfers or measure weight-bearing during gait training. Pair via Bluetooth with a phone to track progress using Rifton's Gait Tracker app.

**Note:** The scale must be included in the original purchase. It cannot be added later.



## Odometer

No more counting tiles! The odometer easily records walking distance in feet or meters.



## Front handle

Allows caregivers to maneuver and guide the E-Pacer, with up/down switch for body support height adjustment.



## Switch pole

Lets caregivers operate the electric lift without bending down. Reversible design gives the client control when desired.



## Arm prompts

Height, rotation, angle, and positioning adjustments.



## Arm platforms

A simple arm prompt alternative for positioning.

**Padded surface dimensions:** 13½" x 5½"



## Hip positioner

Encourages proper positioning and natural movement during gait training. Provides abduction support with optional comfort pad.



## Pelvic support

Available in three sizes, the pelvic support assists with sit-to-stand transfers, gait training, and fall prevention—ideal for early mobility recovery.



## Thigh straps

For seated transfers—machine washable and available in narrow or wide styles.



### Ankle prompts

Prevent scissoring and control stride length with padded straps and spring adjusters.



### Accessories tote

Keeps all E-Pacer components together. Large, durable, and attractive design.

## Need support?

The Rifton website offers video tutorials, training resources, sample letters of medical necessity, evidence-based articles, and more.




Scan here to order, request a quote, or see detailed product information.

## Frame color options



Red

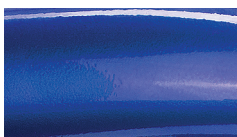

Blue



Black



Lime



Pink

