



Tricycle

Adaptive tricycles give non-ambulatory people the joy and freedom of independent movement. **The Rifton Tricycle** combines four decades of expertise with therapists' insights and cutting-edge technology to create the ultimate ride. It's sleek, rugged, adaptable, and a whole lot of fun.



Scan here to order, request a quote, or see detailed product information.

- Accessible and functional design
- Tool-free adjustments
- Therapeutic benefits including improved leg strength, posture, balance, and range of motion



Key features

1 Seat adjustability

The seat glides smoothly up and down on a low-friction rail with height markings in increments.



2 Durable materials

Built to last with premium materials: carbon fiber drive belts, high traction solid urethane tires, and automotive-grade powder-coated frame finish.



The Tricycle in use



3 Variable resistance and parking brake

The adjustable brake can provide resistance or lock the wheels for safe transfer.



4 Freewheel option

The freewheel option allows caregivers to disconnect the pedals from the drivetrain so the tricycle can be pushed by the caregiver while the pedals remain stationary.



Key features

5 Backrest adjustments

Position the backrest vertically or at a $7\frac{1}{2}^\circ$ or 15° recline.



6 Easy transport

Remove rear wheels and front assembly, then flip down the backrest for compact transport.



The Tricycle in use



Tricycle options



Scan for sizing and pricing.

Frame size options



Large



Medium




Small

Frame color options


Small and medium Tricycle colors




Blue

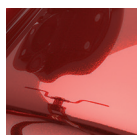

Red

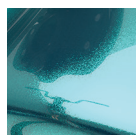

Teal

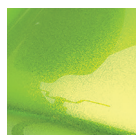

Lime




Pink


Large Tricycle colors



Blue

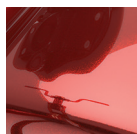

Red

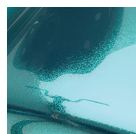

Teal



White



Black

Handlebar options



Handlebars

Choose between stylish conventional handlebars, or loop handlebars for more adjustability.

Backrest options



Trunk support

Stabilizes riders with low trunk control. Includes butterfly harness for front support, with optional lateral supports and headrest.



Backrest pad

For clients who need less support.

Tricycle components



Seats

Soft foam for comfort. Includes a seat belt.

Small: 8" wide x 8" long

Medium: 13" wide x 10" long

Large: 17" wide x 12" long



Headrests

Winged or contoured headrests adjust horizontally and vertically for optimal positioning. Tool-free attachment to trunk support.



Laterals

Slide in trunk support tracks with independent vertical, width, and angle adjustment. Swing out of the way for easy transfers.

Small: 4" x 5"

Large: 5" x 6"



Push handle

Allows caregivers to assist riders and control tricycle speed. Available with or without caregiver brake.



Rear steer

Lets caregivers assist with steering from behind, while giving riders a sense of independence. Available with or without caregiver brake.



Hand anchors

Offer comfortable support for users with limited hand control or grip strength.



Front handbrake

Provides speed control for riders. Adjustable for left or right-hand use. Includes parking lock feature.



Front guide bar

Allows caregiver guidance and enables eye contact between the caregiver and rider.



Abductor

Prevents leg scissoring by separating knees and controlling spasticity. Forward/back adjustment with depth markings.



Abductor with adduction straps

Adduction straps attach to the abductor pad, creating comfortable knee boundaries that maintain proper body and joint alignment.



Pelvic harness

Secures hips and upper thighs without abdominal pressure, providing stable positioning for postural control and proper spinal alignment. Controls extensor spasticity that causes sacral sitting.



Pedal raise kit

Accommodates leg length differences with stackable blocks. Use one or two blocks per side as needed.

Block thickness: ½"



Stationary stand

Converts tricycle into an indoor trainer for pedaling practice and all-weather exercise.



Mounting bar

Clips to the handlebar center and enables device mount attachment. Ideal for communication devices.



Front basket

A convenient storage space for small items.



Outdoor cover

Weather-resistant cover protects your tricycle during outdoor storage.



Freewheel

Allows caregivers to disconnect the pedals from the drivetrain so the tricycle can be pushed by the caregiver while the pedals remain stationary.

Need support?

The Rifton website offers video tutorials, training resources, sample letters of medical necessity, evidence-based articles, and more.



Scan here to order, request a quote, or see detailed product information.